



# CAP THEM ALL!

**Now available!** LVL fully-automated capper / decapper.  
The best choice for high throughput sample storage workflows.

- Available in SBS 96, 48 and 24 formats for all models of LVL SAFE® 2D barcoded tubes
- Automated height detection
- Capping/Decapping time respectively approx. 45 seconds

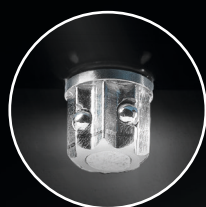

Latest cap driver technology: Cap drivers with spring-loaded metal balls optimally adapt to the caps and permit 100% leak-tight and material-saving sealing of the tubes without tilting.

 Technical specifications located on back

Find out more:  
**LVL cappers online**



**SAFE**   
Store Aliquots For Ever


# Order now.

[lvl-technologies.com/en/devices/capper--decapper](http://lvl-technologies.com/en/devices/capper--decapper)

Innovative and low-maintenance rack-to-head technology:  
It is not the complex driver head moves to the rack,  
but the rack is lifted to the head.

- Can be operated as an integrated and standalone device
- Opens and closes a full rack of tubes at once
- Rack barcode recognition
- Touchscreen panel provides easy navigation through the system
- Caps can be hold or placed on a cap rack
- Ejects caps safely
- Very small footprint 408 x 234 x 388 mm

## Optional accessory LVL Cap Rack

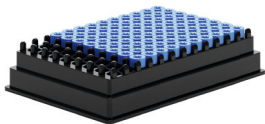


Conveniently store tube caps  
on our cap racks.

SBS 96: RCK-CR96-P-L  
SBS 48: RCK-CR48-P-L  
SBS 24: RCK-CR24-P-L



**96XT/96IT**  
fully-automated capper

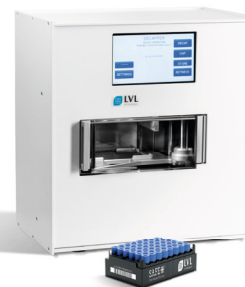

**48XT**  
fully-automated capper



**24XT**  
fully-automated capper

<b>Rack format</b>	<b>SAFE® SBS 96</b>	<b>SAFE® SBS 48</b>	<b>SAFE® SBS 24</b>
<b>Article number</b>	External thread: <b>CDC-96CH-XT-DD</b> Internal thread: <b>CDC-96CH-IT-DD</b>	<b>CDC-48CH-XT-DD</b>	<b>CDC-24CH-XT-DD</b>